
Twin V Low Profile LSO

Warnings and Instructions: Review carefully, proper application is required

⚠**Warning:** This device will not prevent or eliminate risk of injury. Do Not Overtighten. If swelling, pain, skin irritation, or an unusual reaction occurs, discontinue use immediately and consult your medical professional.


Care: Hand wash using mild soap. Rinse thoroughly. Air dry only. Do not tumble dry.

1.

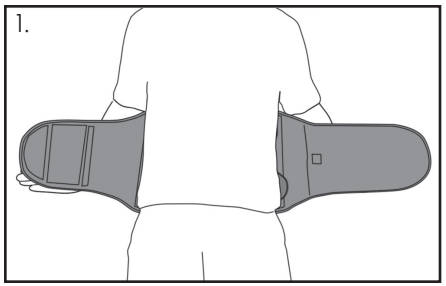


1. Release side pull straps, and position the brace with the back posterior panel centered over the lumbar area.

2.

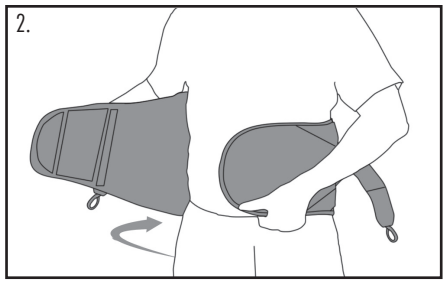


2. Stretch out right side of brace to its furthest extension. Snugly secure right panel over left panel.

3.

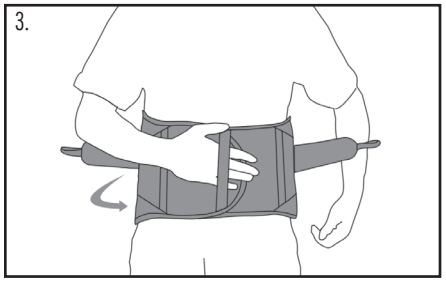


3. Center brace. The front plastic panel should be centered over the abdomen.

4.

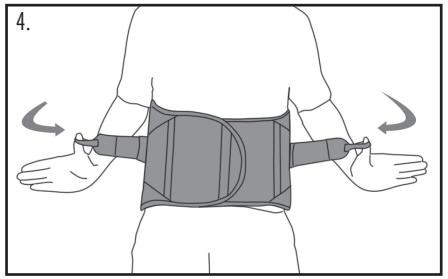


4. Tighten the pull straps and secure.

5.

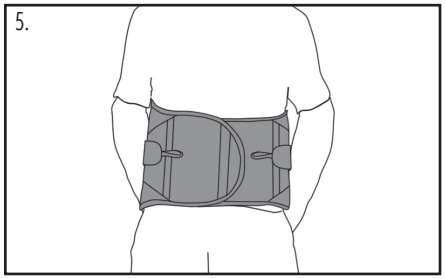

5. Finished application.

6.

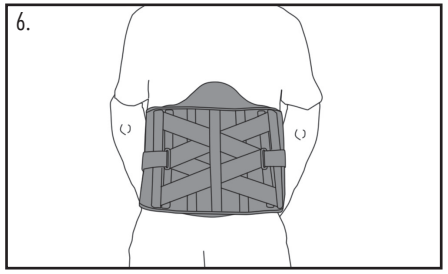

6. Finished application back.

View Application Video

HELY&WEBER
2017 Hely & Weber

409


Twin V Low Profile LSO

Warnings and Instructions: Review carefully, proper application is required

⚠Warning: This device will not prevent or eliminate risk of injury. Do Not Overtighten. If swelling, pain, skin irritation, or an unusual reaction occurs, discontinue use immediately and consult your medical professional.


Care: Hand wash using mild soap. Rinse thoroughly. Air dry only. Do not tumble dry.

1.

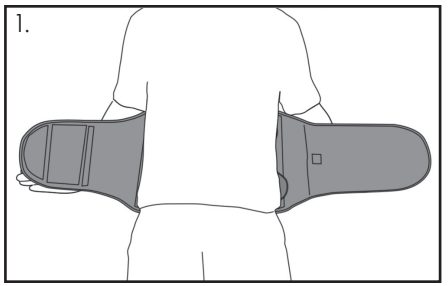


1. Release side pull straps, and position the brace with the back posterior panel centered over the lumbar area.

2.

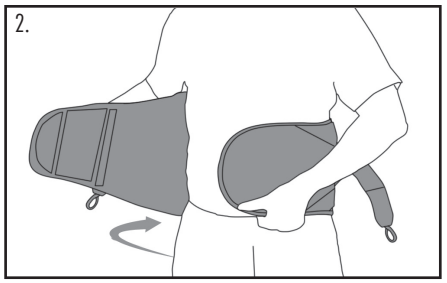


2. Stretch out right side of brace to its furthest extension. Snugly secure right panel over left panel.

3.

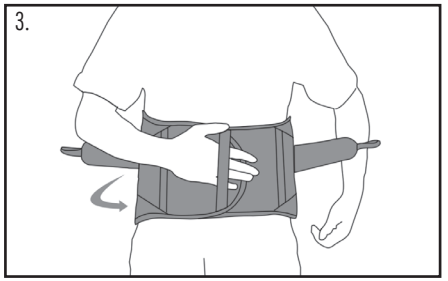


3. Center brace. The front plastic panel should be centered over the abdomen.

4.

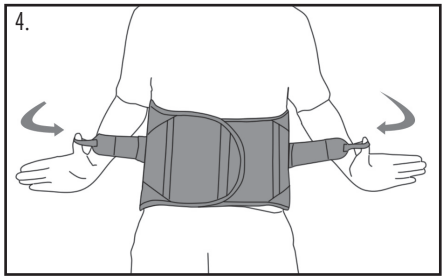


4. Tighten the pull straps and secure.

5.

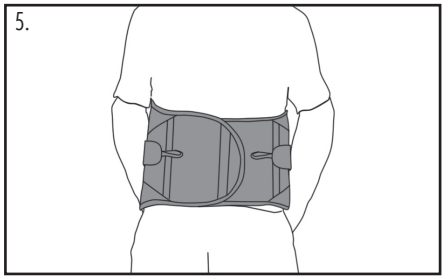

5. Finished application.

6.

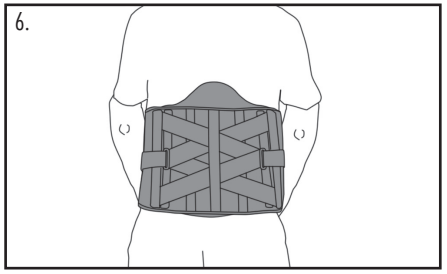

6. Finished application back.

View Application Video

HELY&WEBER
2017 Hely & Weber

Rev D 10/26