



Warnings and Instructions: Review carefully, proper application is required


⚠ Warning: This device will not prevent or eliminate risk of injury. Do Not Overtighten. If swelling, pain, skin irritation, or an unusual reaction occurs, discontinue use immediately and consult your medical professional. This device should not be worn by persons with known allergies to neoprene. Care: Hand wash using mild soap. Rinse thoroughly. Air dry only. Do not tumble dry.

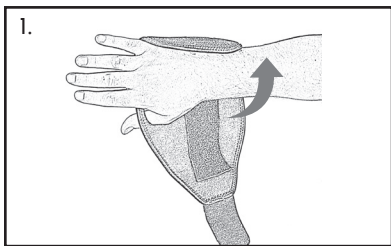


1. Place thumb in opening, pull and fasten inside strap.

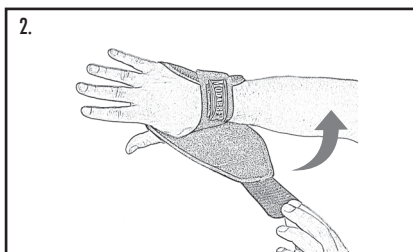


2. Pull and fasten outside strap.

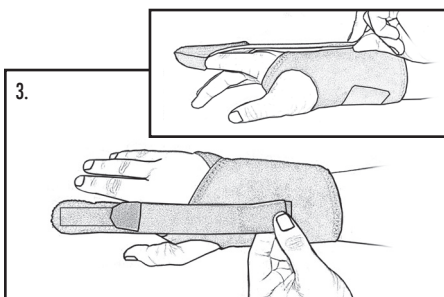


3. Place fully extended finger in finger strap, pull and fasten strap to brace.
DO NOT OVER-EXTEND THE FINGER

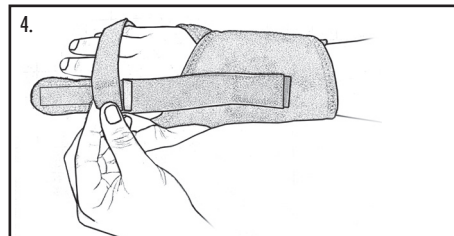

4. Attach small strap to top of finger strap near the hand.

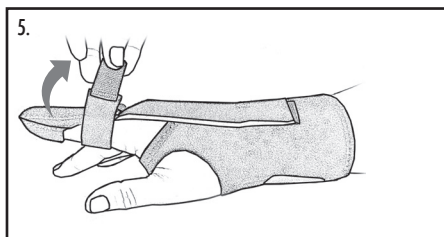

5. Pull small strap around finger and fasten.
DO NOT OVER-TIGHTEN

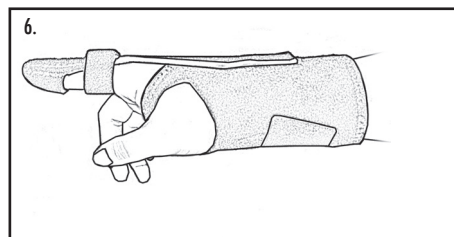


6. Finished application.


U.S. Patent No. 6142966, 7135006

Warnings and Instructions: Review carefully, proper application is required


⚠ Warning: This device will not prevent or eliminate risk of injury. Do Not Overtighten. If swelling, pain, skin irritation, or an unusual reaction occurs, discontinue use immediately and consult your medical professional. This device should not be worn by persons with known allergies to neoprene. Care: Hand wash using mild soap. Rinse thoroughly. Air dry only. Do not tumble dry.

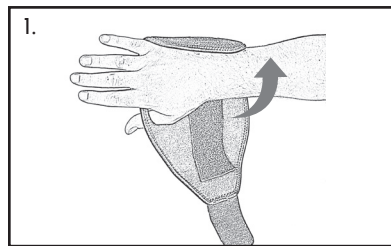


1. Place thumb in opening, pull and fasten inside strap.

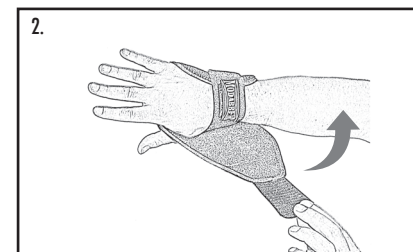


2. Pull and fasten outside strap.

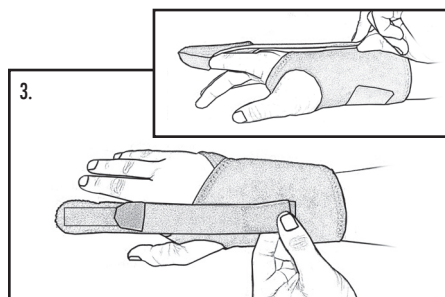


3. Place fully extended finger in finger strap, pull and fasten strap to brace.
DO NOT OVER-EXTEND THE FINGER

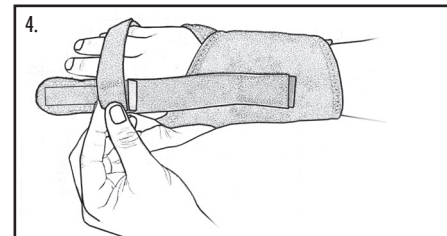

4. Attach small strap to top of finger strap near the hand.

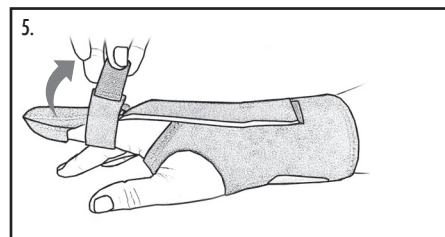

5. Pull small strap around finger and fasten.
DO NOT OVER-TIGHTEN

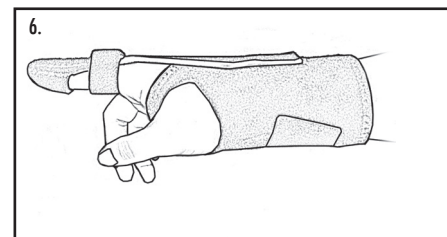

6. Finished application.

U.S. Patent No. 6142966, 7135006